

# R3

 | Knee  
Massager

Instruction Manual



Vibration  
Massage



Hot  
Compress



Infrared  
Therapy



2-Hour  
Battery Life



**ELEEEELS**


## R3 at a Glance



- 1) User-friendly and smart touch panel control
- 2) Built-in 11 red lights
- 3) Infrared laser lamp
- 4) Multi-frequency vibration massage
- 5) Velcro straps
- 6) USB type-c charging port

### Specifications

- Rated Voltage: 3.7V
- Input: DC 5V/2A
- Battery Capacity: 3000 mAh
- Rated Power: 10W
- Weight: 720 grams

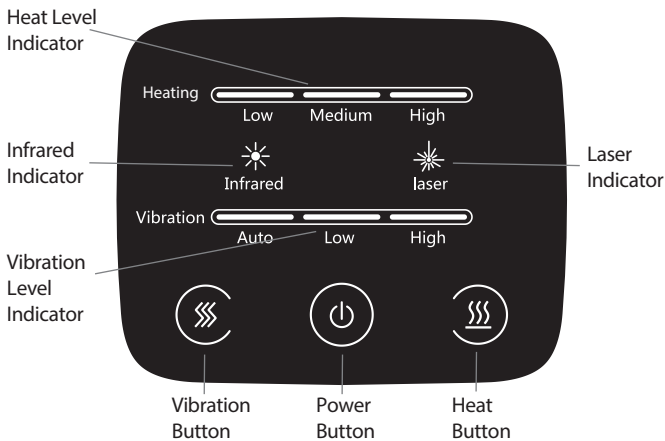

## **INSTRUCTIONS**

PLEASE READ PRIOR TO USE

### **CHARGING**

- 1) Connect USB Type-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use an adapter supporting Type C (5V/2A) to fully charge in approximately 3 hours.
- 3) When the power indicator flashes, please recharge the battery.

### **OPERATING**

- First, use the strap to hold the device on the knee joint, and adjust the tightness.
- Then, press the power button for 2 seconds to turn on the device. The default setting is Auto Mode (auto vibration, low heat level, infrared, laser).
- Press the vibration button repeatedly to adjust the vibration level during operation.
- Press the heat button repeatedly to adjust the heat level during operation.
- Press the power button repeatedly to turn on and turn off infrared and laser during operation.
- Press the power button for 2 seconds to turn off the device.
- When the device stops operation after a 10-minute session, press the power button to start a new one.

### **MAINTENANCE, CLEANING & STORAGE**

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB Type-C cable supplied by the ELEELS for recharging.

### **WARNINGS & CAUTION**

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

## **SAFETY INSTRUCTIONS**

**DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:**

- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

## **DISCLAIMER**

Please do not use the ELEEELS R3, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

## **WARRANTY**

**IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE**

**THIS WARRANTY EXCLUDES:**

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEELS instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:  
<http://www.eleeeels.com/im>