

R6

Eye Care Massager



ELEEELS

Instruction Manual



Soft and
Skin-Friendly
Nano Silkin™
Lining



Graphene
Heating



Air
Acupressure



Vibration
Massage



180°
Folding


R6 at a Glance



- 1) Graphene Heating
- 2) Air Acupressure
- 3) Vibration Massage
- 4) Adjustable elastic strap
- 5) USB-C Charging Port
- 6) USB-C Cable

Specifications

- Rated Voltage: 3.7V
- Input: DC 5V/1A
- Battery Capacity: 1200 mAh
- Rated Power: 5W
- Weight: 318 grams
- Operating Time: up to 2 hours of use per charge

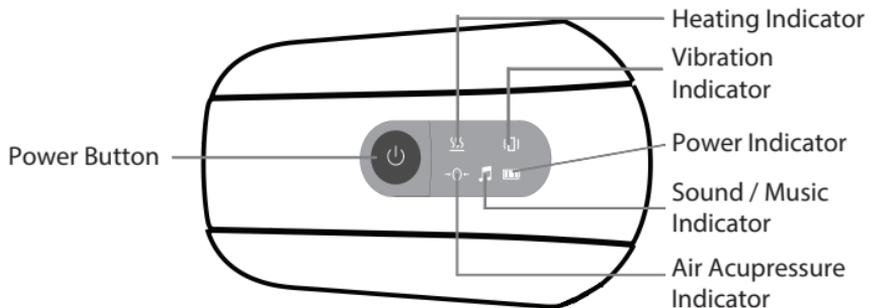

INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

- 1) Connect USB-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use an USB adapter supporting 5V/1A to fully charge in approximately 2 hours.
- 3) During charging, the power indicator flashes slowly. When the power indicator stops flashing, the device is fully charged.

OPERATING - BUILT-IN NATURAL SOUND

- Close your eyes.
- Put on the device. Adjust the strap to the most comfortable position.
- Press the power button for 2 seconds to turn on the device. The default setting is Vitality Mode with built-in natural sound.
- Press the power button repeatedly to switch between 4 modes.
 - **Vitality Mode: Acupressure + Heating + Vibration + Natural Sound**
 - **Rejuvenation Mode: Acupressure + Heating + Natural Sound**
 - **Soothing Mode: Acupressure + Heating**
 - **Sleep Mode: Acupressure + Natural Sound**
- Please note that natural sound is off in Soothing Mode.
- **Press the power button for 3 times consecutively to turn off the natural sound. Press the power button for 3 times consecutively again to resume the natural sound.**
- When the device stops operation after a 10-minute session, press the power button for 2 seconds to start a new one.
- Press the power button for 2 seconds to turn off the device.
- After massage, keep your eyes closed for a few seconds before opening them slowly.
- When the battery level is lower than 10%, the power indicator flashes. Please recharge the battery.

OPERATING - SMARTPHONE MUSIC VIA BLUETOOTH CONNECTION

- To connect R6 with your smartphone, go to the Bluetooth setting of your smartphone, and select "ELEEELS R6 Eye Massager". The sound / music indicator turns blue when pairing is successful, and you can listen to music from your smartphone.
- When R6 connects with your smartphone, smartphone music is on for all 4 modes.
- Press the power button for 3 times consecutively to turn off the smartphone music. Press the power button for 3 times consecutively again to resume the smartphone music.

MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB-C cable supplied by the ELEEELS for recharging.

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.

- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before use.
- Do not use the device while driving.
- Remove the contact lens before using the device.
- Clean and dry your eye area. Make sure it is not oily or sweaty before using the device.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

DISCLAIMER

Please do not use the ELEEELS R6, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: retinal damage, recent surgery, epilepsy, glaucoma, cataract, apoplexy, eye injury, eye inflammation, retinal detachment, pink eye, sensitive skin, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEELS instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:
<http://www.eleels.com/im>